

Evaluation of Prolevo Heelsafe Pressure Relief Pad

DUNCAN BAIN CONSULTING, 22 Gypsy Lane, Kings Langley, Herts, WD4 8PR Tel: 01923 290033

Introduction

This report evaluates the performance of the Prolevo HeelSafe Pressure Relief Pad, which is a leg-section inflatable overlay aimed at reducing interface pressures on the heels. High interface pressures are known to be associated with the onset of pressure ulcers, and these commonly occur on the heels. The report also evaluates the effect of the product on the skin microclimate, as this is also known to be a factor in the aetiology of pressure ulcers.

Methods

Peak interface pressures

These are assessed using the UCL Phantom (developed by the RAFT Institute), a full technical description of which is published in the scientific literature.¹ This is a life-sized articulated dummy with soft tissues, and bony prominences within. The Phantom has an automated positioning system, which places it in exactly the same way on every mattress. Pressure measurements are made using a highly flexible pressure-mapping array, to locate the peak pressures (which occur in different anatomical regions on different mattresses). The surface of the Phantom is warmed to 35°C using special heated and temperature-controlled skin.

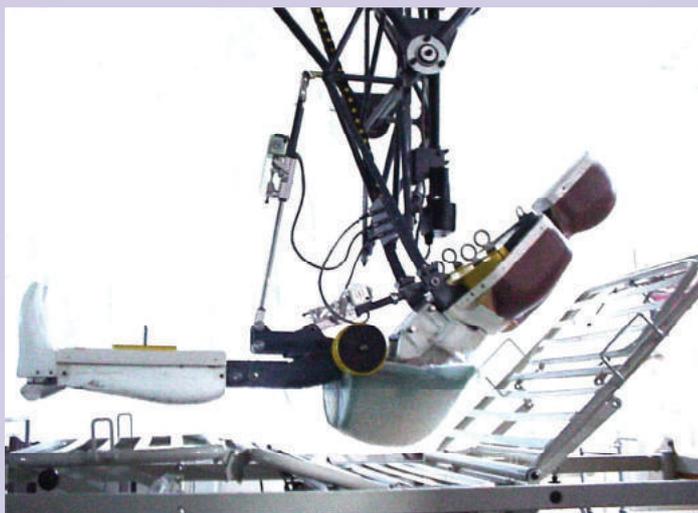


Figure 1 UCL Phantom

Multiple measurements are made, to obtain confidence intervals for the peak pressures in the pelvic and heel regions. Low peak interface pressure is deemed to be the most valid measure of pressure reducing properties according to current evidence at the time of publication.²

Aims

- To evaluate the effect of the HeelSafe overlay on peak interface pressures on the heel and lateral malleolus.
- To establish the optimum pressure setting for the HeelSafe product at a range of different leg weights representing different patient bodyforms.
- To evaluate the heat and water vapour transfer properties of the HeelSafe overlay.

Pressure maps reveal visually much information besides peak pressure about the way pressure is distributed. A picture of the pressure map is therefore also provided, to allow readers to judge features of the pressure distribution that may be of particular interest to them (eg contact area). The pressure maps are provided with a sidebar scale to the colour map.

In the case of the HeelSafe overlay, the interface pressure under the heels is the relevant measurement. The phantom has a calcaneus insert based on a 50th percentile female heel, with a thin covering of gel material based on rheological measurements of hospital patients.

After positioning on the bed, the knee joint and ankle joint is released to allow natural movement.

Pressures were compared on a Softform Premier mattress, with and without the HeelSafe overlay.

Interface pressures were measured using three different setting valves (10mmHg, 20mmHg, and 30mmHg inflation settings), at different body weights of 40kg, 60kg, 80kg, 100kg, 150kg, and 200kg.

Since a side-lying position is frequently adopted, this position was modelled using a live female subject age 71, weighing 62kg, and of height 1.67 m.

Results
Heel pressures on flat bed with HeelSafe

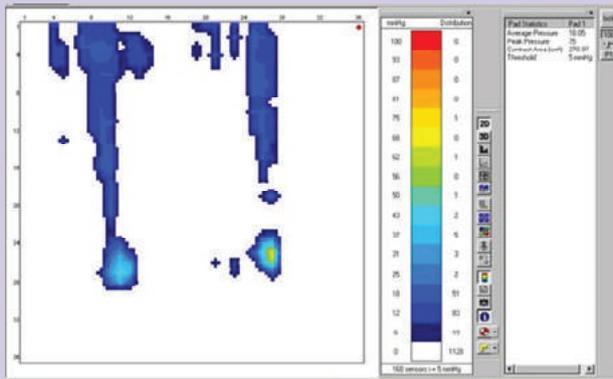


Figure 3 Interface pressure on heels, 20mmHg inflation, 60kg body type

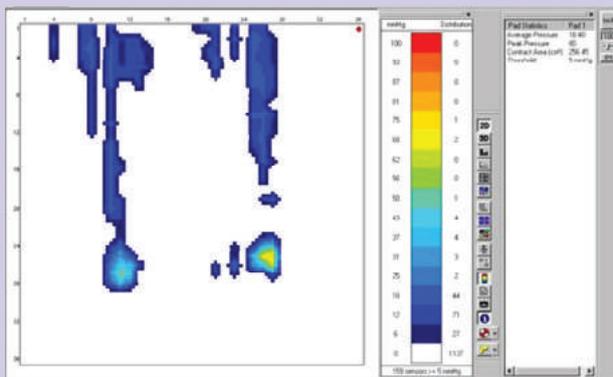


Figure 4 Interface pressure on heels,with HeelSafe , 30mmHg inflation, 60kg body type

Leg break activated with HeelSafe

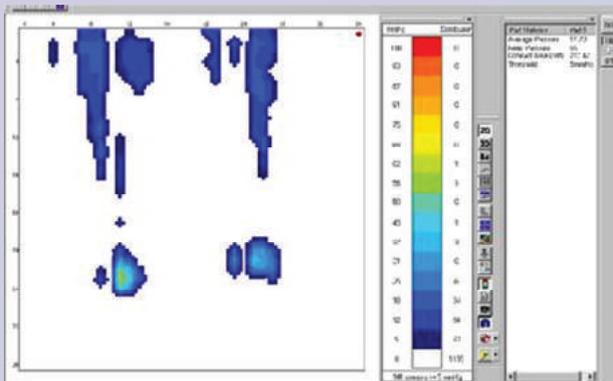


Figure 5 Interface pressure on heels,with HeelSafe , 10mmHg inflation, 60kg body type with leg break


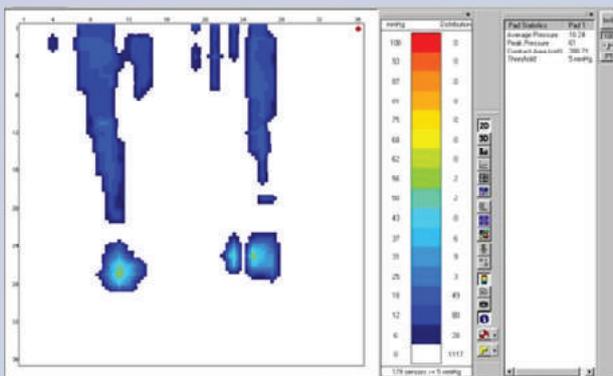


Figure 7 Interface pressure on heels,with HeelSafe , 30mmHg inflation, 60kg body type with leg break

Figure 2 Interface pressure on heels,with HeelSafe , 10mmHg inflation, 60kg body type

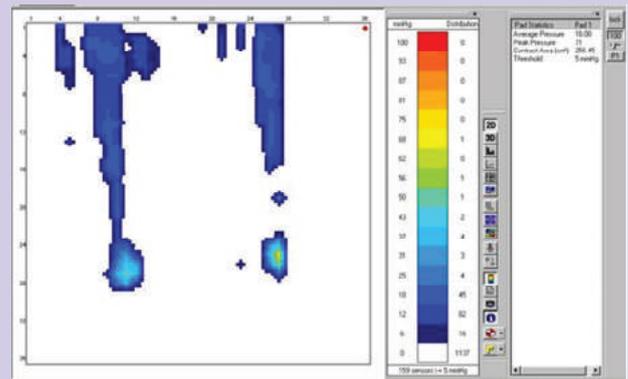

Example pressure maps are shown in figures 2-4 of the effects of inflation pressure on the pressure distribution with Heelsafe. The large array of variables (6 different patient weights x 2 different knee break angles x 3 different inflation pressures x 3 different surface combinations x 2 anatomical sites) gives too many permutations to show every combination tested. The pressure maps shown are selected for illustration. Overall results are collated graphically. It can be seen that increasing the inflation pressure reduces the overall contact area, as expected. However, this does not result in an increase in peak interface pressure. In this orientation, we see that the middle pressure value gives the lowest value of interface pressure. No 'bottoming out' was observed with the lowest pressure, so it is assumed that the slightly elevated pressure was a result of increased surface 'hammock' tension associated with greater penetration.

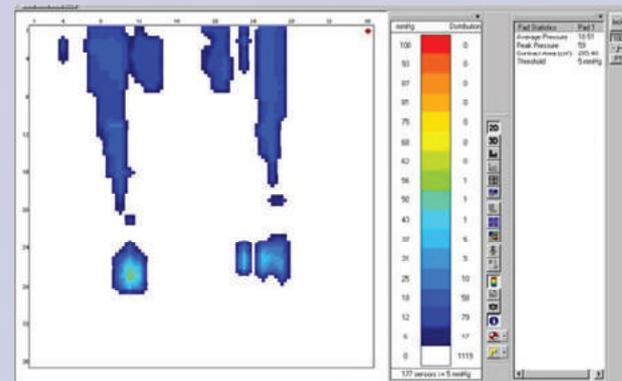


Figure 6 Interface pressure on heels,with HeelSafe , 20mmHg inflation, 60kg body type with leg break

Figures 5-7 show the same procedure with the bed adjusted to allow a 10 degrees flexion of the knee. Again, the middle pressure setting proves to be the optimum for pressure reduction. Generally, lower interface pressures were observed with the bed adjusted to allow knee flexion. This is as expected, since the HeelSafe product lifts the foot slightly from the level of the mattress, causing more leg weight to be transferred through the knee in extension. Flexion at the knee removes moment transferred through the joint, thus removing this extra weight.

Figures 8 and 9 show the peak pressures for all body weight models. It can be seen that, although a trend exists for higher peak pressures at different weights, the middle pressure setting is optimal in most cases. At very light weight, 40kg, the lowest pressure setting gave a lower peak interface pressure, but this was marginal.

It is recommended that the md setting is used across all patient weights, to simplify the decision making process.

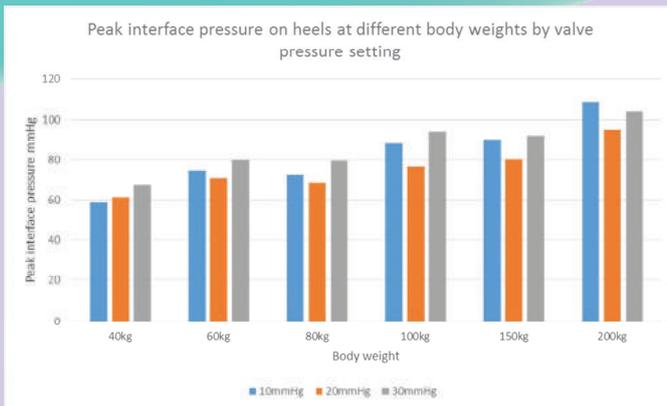


Figure 8 Peak pressures for all body weight models

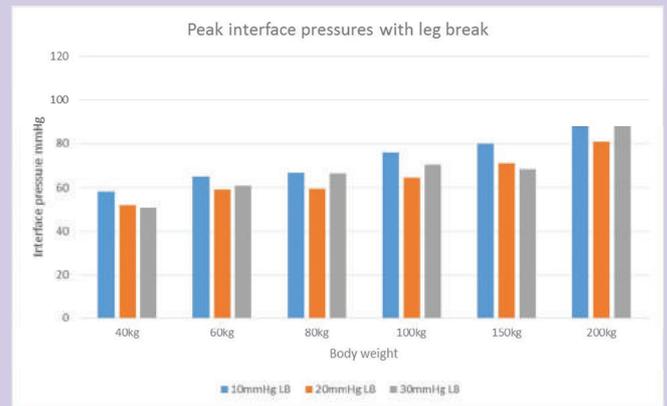


Figure 9 Peak pressures for all body weight models with leg break.

Side Lying

Side lying pressure was mapped using a live subject, as the biomechanics of this position are too complex to model realistically with a phantom. In this case, the subject lay on the bed on her left side in foetal position with the upper (right) leg extended beyond the pressure map. The pressure map shows the interface between the lateral side of the left foot and the bed.

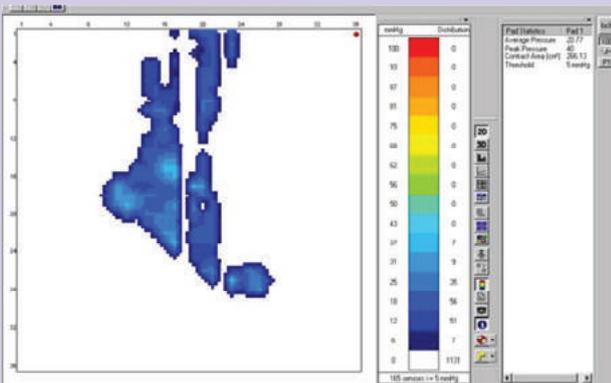


Figure 10 Side lying pressure,with HeelSafe , 10mmHg valve

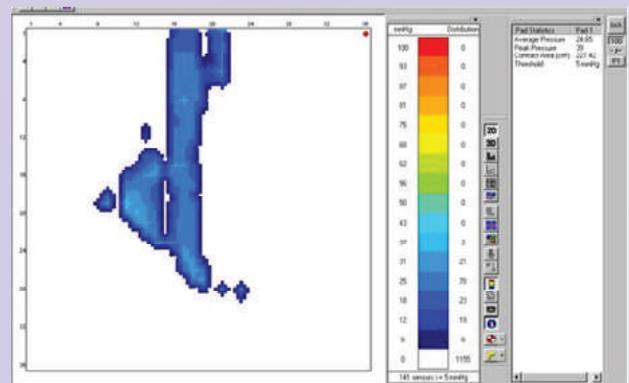


Figure 11 Side lying pressure,with HeelSafe , 20mmHg valve

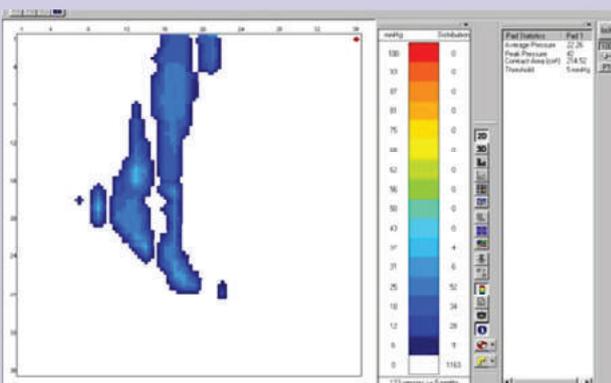


Figure 12 Side lying pressure,with HeelSafe , 30mmHg valve

Also in the side-lying position, the mid-range valve appears to be the optimum setting. The interface pressure seems to be remarkably insensitive to inflation pressure, although it can be seen that there is greater contact area at lower pressure.

Surface comparison: Softform Premier, Domestic mattress (Sealy Pearl Contour), Softform with HeelSafe

Figures 13-15 show the 60kg articulated phantom respectively on the domestic mattress, the Softform Premier, and the Softform Premier with HeelSafe in the optimal setting of 20mmHg inflation. It can be seen that the lowest peak pressures occur on the HeelSafe. This appears to be due to increased contact area, with an almost continuous contact along the Achilles tendon area.

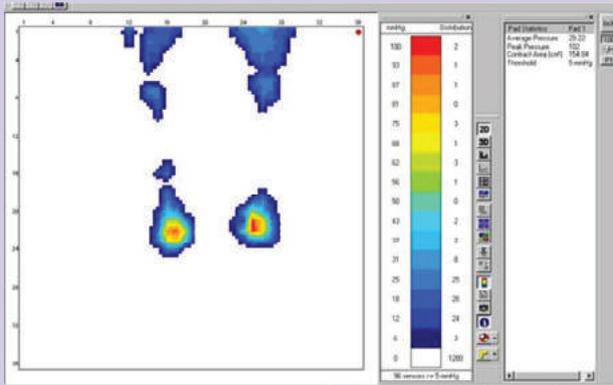


Figure 13 Peak interface pressure, 60kg phantom domestic mattress

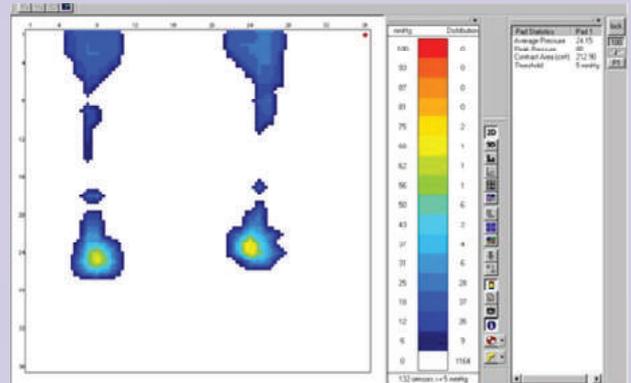


Figure 14 Peak interface pressure, 60kg phantom Softform Premier mattress

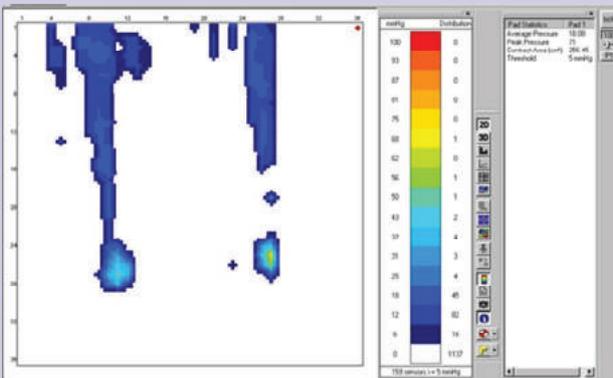


Figure 15 Peak interface pressure, 60kg phantom HeelSafe


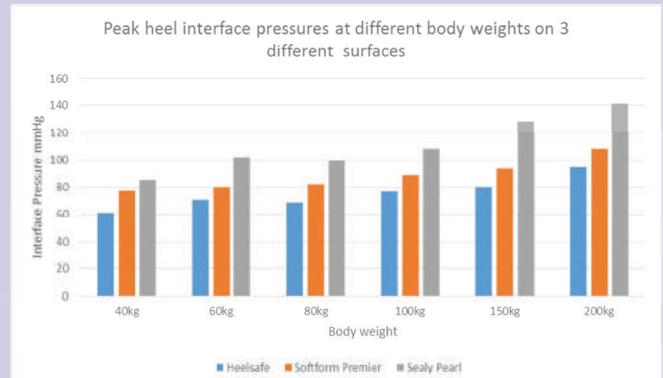


Figure 16 Figure 16 summarises the results at different modelled body weights. It can be seen that in every case the peak pressure was lowest on the HeelSafe compared to the other 2 surfaces.

Microclimate

The ability of a surface to dissipate body heat and moisture makes an important contribution to comfort. Excessively moist conditions at the skin/mattress interface are also known to macerate the skin, exacerbating the risk of mechanical damage to the skin.

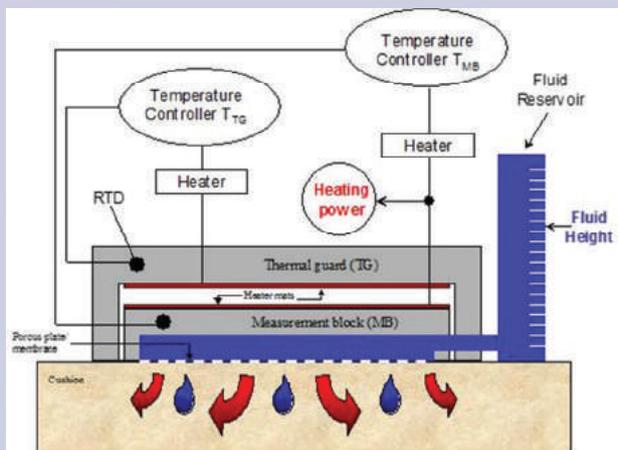

Heat and water vapour transfer test rig schematic.

Tests were done in a controlled environment testing facility with a thermalguarded sweating hot-plate 3. This permits accurate measurements to be made of both heat transfer rates and water-vapour transfer rates through the product.

The hot-plate is maintained at constant temperature and humidity at the interface to the mattress, and losses of heat and water vapour into the mattress are electronically monitored simultaneously.

Test Heat/Vapour permeability with or without cover.

Conclusions

HeelSafe provides effective reduction in heel pressure across all body weights tested, relative to both a Softform Premier mattress and a Sealy Pearl Contour mattress.

Heat and water vapour transfer rates through the product are within expected norms, and no adverse effect on skin microclimate are foreseen, compared to existing pressure reducing mattresses.

Optimum inflation pressure appears to be the mid setting (20mm Hg) of the 3 valves provided, and this seems to be the case across all body weights. This favours a simple approach to inflation, whereby patient weight need not be considered prior to inflating the product.